

A guide to use the flowPeaks: a flow cytometry data clustering algorithm via K -means and density peak finding

Yongchao Ge

May 19, 2021

yongchao.ge@gmail.com

1 Licensing

Under the Artistic License, you are free to use and redistribute this software. However, we ask you to cite the following paper if you use this software for publication.

- Ge Y. et al, flowPeaks: a fast unsupervised clustering for flow cytometry data via K -means and density peak finding, 2012, Bioinformatics 8(15):2052-8.

2 Overview

We combine the ideas from the finite mixture model and histogram spatial exploration together to find the clustering of the flow cytometry data. This new algorithm, in which we called **flowPeaks**, can be applied to high dimensional data and identify irregular shape clusters. The algorithm first uses K -means algorithm with a large K to partition the population into many compact clusters. These partitioned data allow us to generate a smoothed density function. All local peaks are exhaustively found by exploring the density function and the cells are clustered by the associated local peak. The algorithm flowPeaks is automatic, fast and reliable and robust to the cluster shape and outliers. Details can be seen in the paper by Ge (2012).

3 Installation

3.1 All users

When you are reading this and find the R package is already available in the Bioconductor package repository. You can install it directly from Bioconductor.

- **Windows users:** select the menu “Packages” and then click “Select repositories...” and choose “BioC software”. And then select the menu ‘Packages’ click “install R package(s)...” and then look for the package **flowPeaks**.
- **Linux users:** This also works for Windows users. Type the following after you have invoked R

```
> if (!requireNamespace("BiocManager", quietly=TRUE))  
  > install.packages("BiocManager")  
> BiocManager::install("flowPeaks")
```

If this succeeds, congratulations, you can ignore the rest of this section and skip to Section 4.

3.2 Windows Users

Please read section 3.1 to install the R package from Bioconductor before proceeding this. If you have the prebuilt binary of **flowPeaks** zip file, you can install the package by selecting the menu “Packages”, and then “Install packages from a local zip file”, and then point to prebuilt binary of **flowPeaks** zip file.

To build **flowPeaks** from the source by using Rtools is very not straightforward. R novices are not encouraged to try this. Experienced R users need to carefully follow the instruction of the Rtools (<http://www.murdoch-sutherland.com/Rtools/>) and <http://cran.r-project.org/doc/manuals/R-admin.html#The-Windows-toolset>. The GSL library needs to be downloaded from the file **local215.zip** at <http://www.stats.ox.ac.uk/pub/Rtools/goodies/multilib/>. The top folder (**top_path_local**) of the extracted file **local215.zip** should contain three subfolders: **include**, **share** and **lib**. The next step is to modify the file **flowPeaks/src/Makevar.win** as below.

- **PKG_LIBS** += -L(top_path_local)/lib/\$(R_ARCH)/ -lgsl -lgslcblas -lm
- **PKG_CXXFLAGS** += -I(top_path_local)/include

The users are not encouraged to compile their own gsl library by MinGW or Visual Studio. Most likely their own version of gsl library is not going to work.

3.3 Linux Users

To build the **flowPeaks** package from the source, make sure that the following is present on your system:

- C++ compiler
- GNU Scientific Library (GSL)

A C++ compiler is needed to build the package as the core function is coded in C++. GSL can be downloaded directly from <http://www.gnu.org/software/gsl/> and follow its instructions to install the GSL from the source code. Alternatively, GSL can also be installed from your linux specific package manager (for example, Synaptic Package Manager for Ubuntu system). Other than the GSL binary library, please make sure the GSL development package is also installed, which includes the header files when building **flowPeaks** package.

Now you are ready to install the package:

```
R CMD INSTALL flowPeaks_x.y.z.tar.gz
```

If GSL is installed at some non-standard location such that it cannot be found when installing **flowPeaks**. You need to do the following

1. Find out the GSL include location (`<path-to-include>`) where the GSL header files are stored in the sub folder `gsl`, and GSL library location (`<path-to-lib>`) where the lib files are stored. If the GSL's `gsl_config` can be run, you can get them easily by `gsl-config --cflags` and `gsl-config --libs`
2. In the file `flowPeaks/src/Makevars`, you may need to change the last two lines as below:
 - `PKG_CXXFLAGS = -I<path-to-include>`
 - `PKG_LIBS = -L<path-to-lib>`

4 Examples

To illustrate how to use the core functions of this package, we use a barcode flow cytometry data, which can be accessed by the command `data(barcode)`. The barcode data is just a simple data matrix of 180,000 rows and 3 columns. The clustering analysis is done by using the following commands.

```
> library(flowPeaks)
```

```
> data(barcode)
```

```
> summary(barcode)
```

| Pacific.blue | Alexa | APC |
|--------------|----------------|--------------|
| Min. : 153 | Min. : 81.0 | Min. : 92 |
| 1st Qu.:1198 | 1st Qu.: 602.0 | 1st Qu.: 503 |
| Median :1900 | Median : 699.0 | Median :1817 |
| Mean :1817 | Mean : 688.7 | Mean :1565 |
| 3rd Qu.:2522 | 3rd Qu.: 789.0 | 3rd Qu.:2393 |
| Max. :3503 | Max. :1119.0 | Max. :3369 |

```
> plot(fp,idx=c(1,3))
```

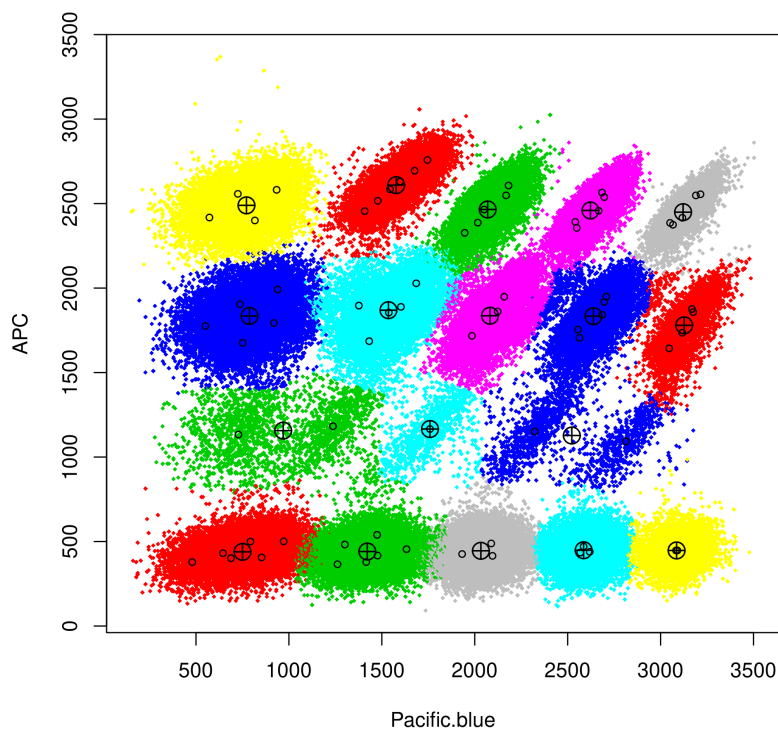


Figure 1: The scatter plot of clustering results for the first column and the third column. Two clusters may share the same color as the automatic color specification is not unique for all clusters. The users may have option to change the col option in the plot function to have their own taste of color specifications. The flowPeaks cluster are drawn in different colors with the centers (\oplus), the underlying K -means cluster centers are indicated by \circ .

```
> fp<-flowPeaks(barcode)
```

```
step 0, set the intial seeds, tot.wss=5.47711e+09
step 1, do the rough EM, tot.wss=3.52259e+09 at 0.652075 sec
step 2, do the fine transfer of Hartigan-Wong Algorithm
tot.wss=3.52028e+09 at 0.951423 sec
```

If we want to visualize the results, we can draw scatter plot for any two columns of the data matrix. The result is shown in Figure 1. Different colors specify different clusters, and the dots represent the center of the initial K -

```
> par(mfrow=c(2,2))
> plot(fp,idx=c(1,2,3))
```

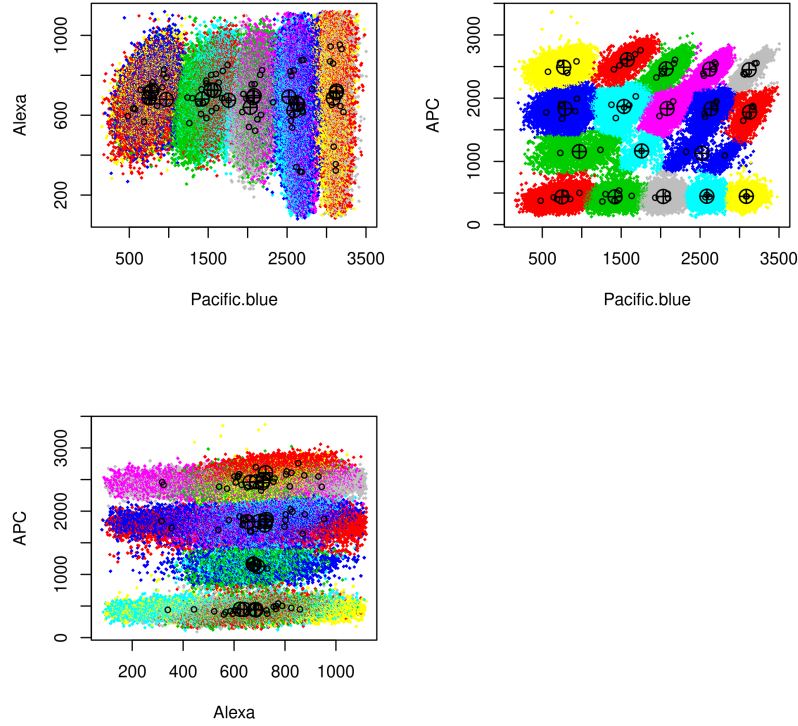


Figure 2: Pairwise plot of all columns

means, and triangles are the peaks found by flowPeaks. We can see all pairwise scatter plots, the R commands and the plot are shown in Figure 2.

We find out that the data displays better clustering in the Pacific.blue and APC. We could choose to redo the clustering just focusing on these two dimensions as shown in Figure 3.

Since the clustering is on two dimension, the plot in Figure 3 is able to give you the boundary of the final clusters (in bold line) and the Voronoi boundary of the K -means. These boundaries are not visible for the data clustering in higher dimension as the projection onto 2D will have many clusters overlapped.

We can evaluate the cluster performance with the gold standard clustering stored in barcode.cid as shown in Figure 3's legend by using the `evalCluster()` function.

We could remove the outliers that are far away from the cluster center or

```

> fp2<-flowPeaks(barcode[,c(1,3)])
> plot(fp2)

> evalCluster(barcode.cid,fp2$peaks.cluster,method="Vmeasure")

[1] 0.9983352

```

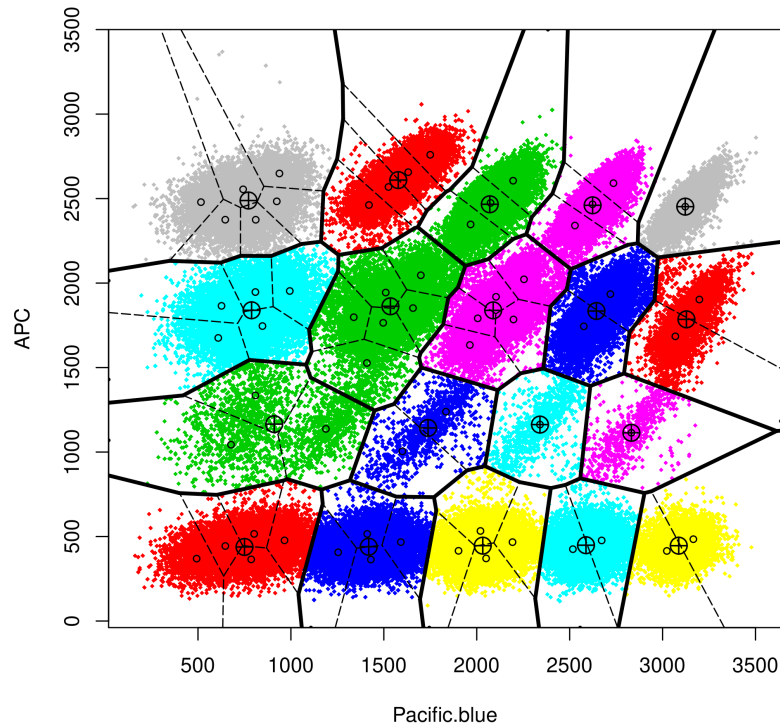


Figure 3: Clustering results based on only the Pacific.blue and APC

can not be ambiguously assigned to one of the neighboring cluster as shown in Figure 4.

5 Reading the data from flowFrame

If a user has the data that is of class `flowFrame` in the `flowCore` package, he just needs to use the transformed data of the expression slot of the `flowFrame` where only channels of interest are selected. See Figure 5 for the details

```

> fpc<-assign.flowPeaks(fp2,fp2$x)
> plot(fp2,classlab=fpc,drawboundary=FALSE,
+      drawvor=FALSE,drawkmeans=FALSE,drawlab=TRUE)

```

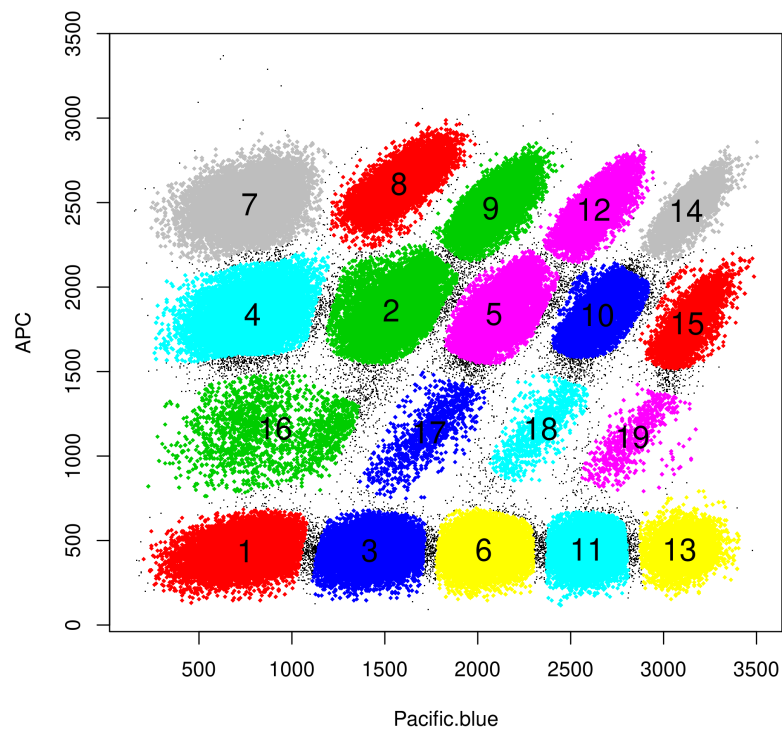


Figure 4: The plot of the barcode after the outliers have been identified as black points

```

> require(flowCore)
> samp <- read.FCS(system.file("extdata", "0877408774.B08",
+ package="flowCore"))
> ##do the clustering based on the asinh transforamtion of
> ##the first two FL channels
> fp.z<-flowPeaks(asinh(samp@exprs[,3:4]))
> plot(fp.z)

```

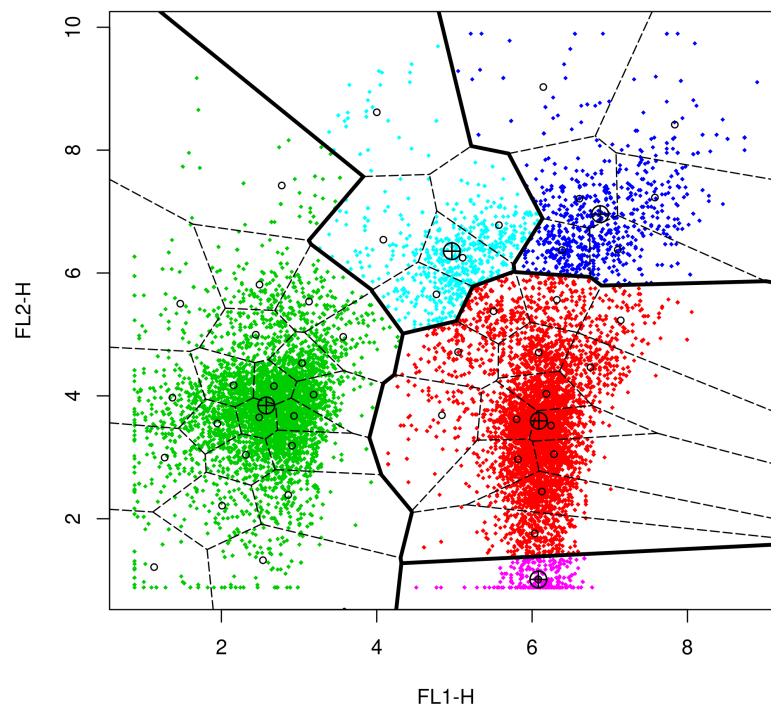


Figure 5: How to use the flowFrame data for the flowPeaks